



DIAGNOSTIC AUTOMATION, INC.

23961 Craftsman Road, Suite D/E/F, Calabasas, CA 91302

Tel: (818) 591-3030 Fax: (818) 591-8383

onestep@rapidtest.com

technicalsupport@rapidtest.com

www.rapidtest.com



See external label



2°C-8°C



Σ=96 tests



Cat # 1740-09

Total Human Lipoprotein a

Lp(a)

ELISA Assay

Cat # 1740-09

Test	Lp(a) ELISA
Method	ELISA: Enzyme Linked Immunosorbent Assay
Principle	ELISA – Solid Phase, Sandwich ELISA
Detection Range	3 µg/dL – 405 µg/dL
Sample	10ul Serum
Specificity	90%
Sensitivity	3 µg /dL
Total Time	120 mins
Shelf Life	12 -18 Months

Intended Use:

To quantify total human lipoprotein A (Lp(a))

Principle of Procedure:

Solid phase capture sandwich ELISA assay using a microwell format.

Shelf Life:

The expiration date for the package and each component is stated on the label(s). Store all components at 2-8°C with the exception of the standard, which should be stored frozen.

Patient and Standard Dilutions:

Dilute each serum or plasma specimen to be tested 1:400 with the Lp(a) specimen diluent provided. (Serum specimens with high Lp(a) levels should be diluted more than 1:400 for accurate Lp(a) determination.)

***Note: A pre-dilution using PBS (phosphate buffer) may be done followed by a final dilution in specimen diluents to bring the serum or plasma final dilution to 1:400.**

Construct a standard curve as follows:

- a) Perform a series of at least four, twofold dilutions of the 1:400 standard. Use the specimen diluent alone as a blank or zero control.
- b) Use the declared value on the vial of Lp(a) standard to calculate the values on the standard curve.

Materials Supplied:

1. Goat Anti-Human Lp(a) coated microwell strips 12x8 with plastic frame
2. Lp(a)N Conjugate – 12mL
3. Lp(a) standard (diluted 1:400) – 1 mL
4. TMB/peroxide substrate color developer –12Ml
5. Lp(a) specimen diluent – 60mL
6. Sulfuric acid termination reagent (0.5N) –12mL
7. 15 X Wash buffer concentrate – 60mL

Limitations of the Procedure:

No single assay should be used as the only basis for arriving at a diagnostic conclusion. For research use only.

Dynamic Range:

3 µg/dL – 405 µg/dL

Reproducibility:


C.V. 4%-8% depending upon the region of the standard curve.

Assay Procedure:

*Allow each reagent to reach room temperature before use

1. Add 100uL of *diluted* specimen or standard to each microwell
2. Incubate at room temperature for 60 minutes
3. Decant and wash each microwell five times with wash buffer (dilute buffer 1:15 with reagent grade water)
4. Add 100uL of anti-human Apo B-100 conjugate to each well
5. Incubate at room temperature for 60 minutes
6. Decant and wash as in step 3
7. Add 100uL of TMB/peroxide substrate and incubate at room temperature for 30 minutes
8. Terminate the reaction with 100uL of 0.5N sulfuric acid
9. Zero the microwell reader at 450nm using the specimen diluent zero control well
10. Determine the optical density (O.D.) of the remaining wells
11. Construct a standard curve using the O.D. values obtained for each of the standards
12. Interpolate the unknowns from the standard curve

Table 1. Lp(a) levels (mg/dl) in centenarians and controls

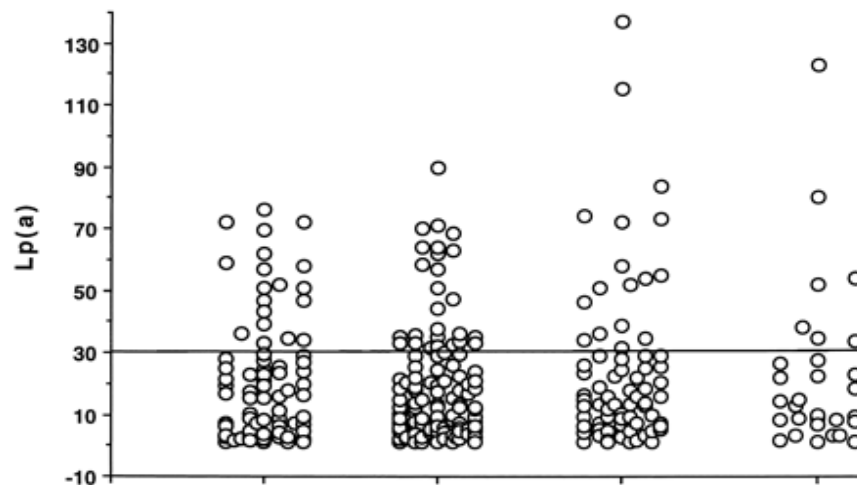

	Centenarians (n=75)	<65 years, randomly selected (n=114)	>65 years, randomly selected	>60 years, healthy selected (n=30)
--	--------------------------------	--	--	--

			(n=73)	
Age range (in years)	100–106	8–64	65–98	61–80
Age mean	100.9 ± 1.4	35.8 ± 11.8	83.5 ± 7.6	71.4 ± 5.5
Lp(a) average	22.4	19.3	23.8	23.0
Lp(a) median	17.2	12.5	15.2	14.2
Lp(a) range	1–76	1–90	1–137	1–123
Log Lp(a)(\bar{x} ±SD)	1.11 ± 0.52	1.06 ± 0.48	1.13 ± 0.51	1.12 ± 0.51
% subjects with Lp(a) >30 mg/dl	25.3	22.8	23.3	23.3
% subjects with Lp(a) <30 mg/dl	74.7	77.2	76.7	76.7

Table is from “Lipoprotein(a) and lipoprotein profile in healthy centenarians: a reappraisal of vascular risk factors” in the FASEB Journal 1998; 12:433-437

Date Adopted	Reference No.
2008-04-28	DA-Lp (a)-2009


DIAGNOSTIC AUTOMATION, INC.
 23961 Craftsman Road, Suite D/E/F, Calabasas, CA 91302
 Tel: (818) 591-3030 Fax: (818) 591-8383
ISO 13485-2003



Revision Date:10/15/2009